



Blood Collection Product Catalog

INNOVATIVE MEDICAL TECHNOLOGIES, INC
OMAHA, NE


866-560-1820
support@InnovativeMedTech.com
www.InnovativeMedTech.com

ULTRAFLO™-L SAFETY BLOOD COLLECTION SET WITH ULTRAFLO™ CANNULA


- Increase 23G flow rate 200%* (60% 21G)
 - Improve patient comfort. Use 23G instead of 21G without compromising flow rate.
 - Proprietary bevel technology requires lower penetration force* for an improved draw experience
- *3rd-party testing compared to leading brand.

ULTRAFLO™-L SAFETY SHIELD


- Same design as the leading brand. No change in technique
- One or two-handed activation

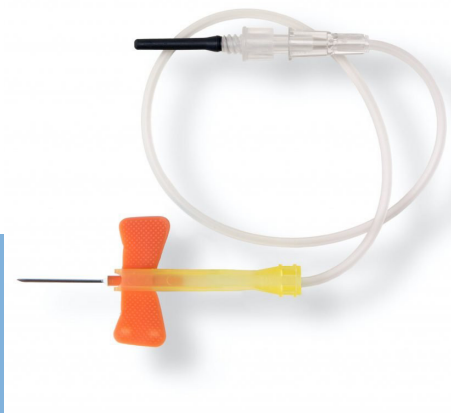

UltraFlo™-L Safety Blood Collection Set 25G. 50/bx. 4bx/cs

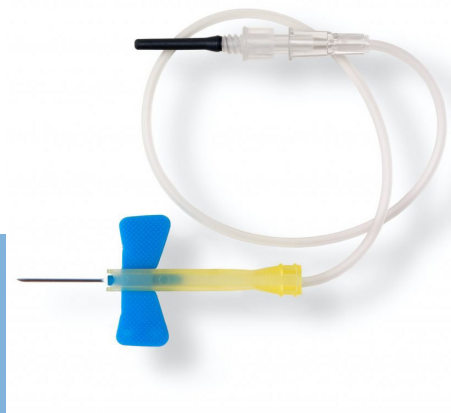

UltraFlo™-L Safety Blood Collection Set 23G. 50/bx. 4bx/cs

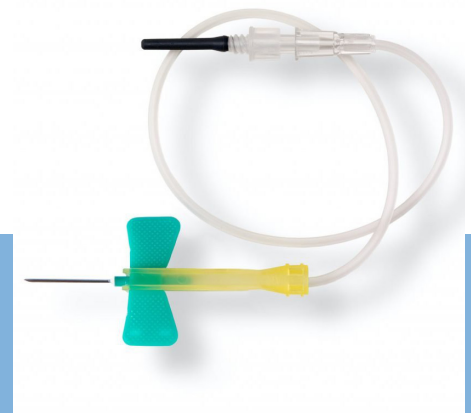

UltraFlo™-L Safety Blood Collection Set 21G. 50/bx. 4bx/cs

Item #	Description
281	UltraFlo™-L Safety Blood Collection Set With MS Luer Adapter, 21 gauge 3/4 in, 12 in Tube.
283	UltraFlo™-L Safety Blood Collection Set With MS Luer Adapter, 23 gauge 3/4 in, 12 in Tube.
285	UltraFlo™-L Safety Blood Collection Set With MS Luer Adapter, 25 gauge 3/4 in, 12 in Tube.
292	UltraFlo™-L Safety Blood Collection Set With MS Luer Adapter, 23 gauge 3/4 in, 7 in Tube.

IDENTICAL SAFETY-SHIELD-ACTIVATION TECHNIQUE

Description	UltraFlo Item #	QTY/BX	BD Item #	QTY/BX
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 21 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	281	50	367281	50
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 23 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	283	50	367283	50
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 25 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	285	50	367285	50
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 23 gauge 3/4 in, 7 in Tube 50/BX 200/CS	292	50	367292	50

DIFFERENT SAFETY-SHIELD-ACTIVATION TECHNIQUE

	COVIDIEN 50/BX	GREINER 50/BX	SMITHS 50/BX	TERUMO 50/BX	KAWASUMI 50/BX	KAWASUMI 50/BX
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 21 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	8881225182	450095	972112	MNSVS 21B30	DMBS- 21G	DMBA- 21G
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 23 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	8881225190	450096	972312	MNSVS 23B30	DMBS- 23G	DMBA- 23G
UltraFlo™-L Safety Blood Collection Set With Luer MS Adapter, 25 gauge 3/4 in, 12 in Tube, 50/BX, 200/CS	8881225208	450099	972512	MNSVS 25B30	DMBS- 25G	DMBA- 25G

UltraFlo™Cannula

- Increase flow rate 200%*
 - Improve patient comfort. Use 23G instead of 21G without compromising flow rate
 - Proprietary bevel technology requires lower penetration force* for an improved draw experience
- *Compared to leading brand

UltraFlo™-L Safety Shield

- Same activation technique as the leading brand --proven and familiar
- One or two-handed activation technique

IMPROVACUTER™ SHATTERPROOF PET TUBES

FDA-approved (k093910) IMPROVACUTER™ Shatterproof PET tubes offer a safe method of blood collection and reduce the potential for tube breakage and specimen spillage, thereby reducing the potential for exposure to bloodborne pathogens.

IMPROVACUTER™ GEL AND CLOT ACTIVATOR TUBES

IMPROVACUTER™ Gel and Clot Activator Tubes offers serum for determinations in chemistry and serology. The tube contains a gel that forms a stable barrier between the serum and the cellular components of the whole blood after centrifugation. They are coated with silicone and micronized silica particles to accelerate clotting. A silicone coating reduces adherence of red cells to tube walls.

The Safety closure has a unique rubber stopper that seals IMPROVACUTER™ blood collection tubes. It is recessed within and covered by a plastic shield to help protect laboratory personnel from contact with blood on the stopper or around the outer rim of the tube. The closure also helps prevent blood from splattering when the tube is opened.



REF#	Material	Tube Size (mm)	Draw Volume (ml)	Closure Type	Cap Color	Label	Additive	Packaging Quantity Box/Case
3050	PET	13 x 100	5.0	Safety	Dark Yellow	Paper	Gel/Clot Activator	100/ 1000
4850	PET	16 x 100	8.0	Safety	Dark Yellow	Paper	Gel/Clot Activator	50/1000

IMPROVACUTER™ is a trademark of Guangzhou Improve Medical Instruments

Improvacuter™ Serum-Gel Blood Collection Tube Cross-Re-

REF #	Material	Tube Size (mm)	Draw Volume (mm)	Closure Type	Cap Color	Additive	Packaging Quantity Box/Case
3050	PET	13x100	5.0	Safety	Gold	Clot Activator/ Separator Gel	100/1000
BD 367986 BD 367989 BD 367977- (4.0 ml)	PET	13x100	5.0	Safety	Gold	Clot Activator/ Separator Gel	
GBO 456018P GBO 456073P	PET	13x100	5.0	Safety	Red/Yellow	Clot Activator/ Separator Gel	
4850	PET	16x100	8.0	Safety	Gold	Clot Activator/ Separator Gel	50/1000
BD 367988 (conventional red/gray stopper)	PET	16x100	8.5	Safety	Gold	Clot Activator/ Separator Gel	
GBO 455071P	PET	16x100	8.0	Safety	Red/Yellow	Clot Activator/ Separator Gel	


PLASTIC BLOOD GAS CAPILLARIES

Our blood gas capillaries are made of fracture-proof, chemically neutral plastics, and, therefore avoid risks of injuries or infections during blood collection. The capillaries are coated with 50 I.U. balanced lithium heparin and can be sealed with closing caps. Special sizes can be manufactured according to your requirements!

HEPARIN SYRINGES FOR VETERINARY

FOR VETERINARY AND LABORATORY USE ONLY:

The Blood Gas Sampling Syringes serve as arterial extraction devices for pH/blood gas, electrolyte, oximetry and metabolite analysis. The devices are coated with approx. 50 I.U. balanced heparin. Sterile, single use. Accessories include clot catchers and air venting caps.

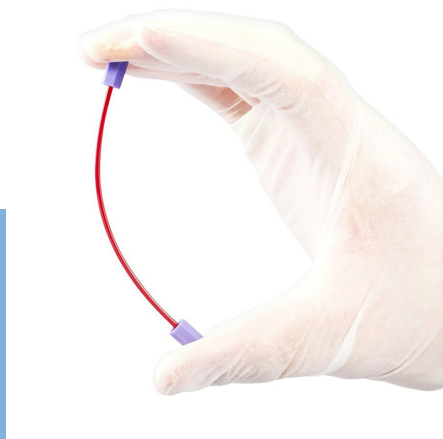

Plastic Blood Gas Capillaries



1 mL Heparin Syringes for Veterinary 50/bx. 20bx/cs

REF	Syringe Description	Box	Case
3999599	2 ml, Luer Slip, Sterile, ~50 I. U. bal Heparin	50 ea	1000 ea
999017	1 ml, Luer Slip, Sterile, ~50 I. U. bal Heparin	50 ea	1000 ea

Plastic Blood Gas Capillaries

Capillary Volume	Length	OuterØ	Pack. Unit ea	Item #
~80 µl	~100 mm	~1.60 mm	100/1000	104110
~90 µl	~80 mm	~1.75 mm	100/1000	101080
~100 µl	~90 mm	~1.75 mm	100/1000	101090
~110 µl	~100 mm	~1.75 mm	100/1000	101100
~130 µl	~120 mm	~1.75 mm	100/1000	101120
~165 µl	~150 mm	~1.75 mm	250/1000	101150-250
~200 µl	~180 mm	~1.75 mm	250/1000	101180-250
~100 µl	~65 mm	~2.35 mm	100/1000	102065
~115 µl	~74 mm	~2.35 mm	100/1000	102070
~140 µl	~90 mm	~2.35 mm	100/1000	102090
~155 µl	~100 mm	~2.35 mm	100/1000	102100
~185 µl	~120 mm	~2.35 mm	100/1000	102120
~230 µl	~150 mm	~2.35 mm	250/1000	102150-250
~80 µl	~54 mm	~1.90 mm	100/1000	102118
~115 µl	~78 mm	~1.90 mm	100/1000	102117
~150 µl	~100 mm	~1.90 mm	100/1000	102116
~170 µl	~100 mm	~2.00 mm	100/1000	105107
Closing Caps		1.75 OD mn	500/1000	110180
Closing Caps		2.0-2.4 OD mn	500/1000	110181
Mixing Wires		15 x 0.7 mn	1000	110009
Mixing Wires		9 x 0.7 mn	1000	110007

ULTRAFLO™-T BLOOD TRANSFER DEVICE

The UltraFlo™-T Blood Transfer Device (Cat. No. 810) is a sterile, pre-assembled and easy-to-use device designed with safety in mind.

It is used for needleless specimen transfer from a syringe to an evacuated tube or blood culture bottle and has a pink color-coded connection to provide easy differentiation from other holder-based products.

ULTRAFLO™-A LUER ACCESS DEVICE

The UltraFlo™-A Access Device (Cat. No. 902) is a sterile, pre-assembled

blood collection tube holder with integrated multi-sample luer adapter. It is compatible with female luer connections or IV ports designed for luer access and has a blue color-coded connection to provide easy differentiation from other holder-based products.



Description	Item #	QTY/CS	BD ITEM#
UltraFlo™-T Blood Transfer Device Female LL (pink) 200/bx. 800/cs	810	800	364810
UltraFlo™-A Luer Access Device Male LL (blue) 200/bx 800/cs	902	800	364902

ULTRAFLO™-V SAFETY BLOOD COLLECTION NEEDLES

The UltraFlo™-V safety multisample blood collection needle features a safety shield that allows for one-handed activation to cover the needle immediately upon withdrawal from the vein and confirms proper activation with an audible click.



UltraFlo™-V Safety Blood Collection Needle 22G 48/bx 10bx/c



UltraFlo™-V Safety Blood Collection Needle 21G 48/bx 10bx/cs

IDENTICAL SAFETY-SHIELD-ACTIVATION TECHNIQUE

Description	UltraFlo™-V Item #	QTY/BX	BD Item #	QTY/BX
UltraFlo™-V Safety Blood Collection Needle 21G 48/bx 10bx/c	607	48	368607	48
UltraFlo™-V Safety Blood Collection Needle 22G 48/bx 10bx/c	608	48	368608	48



INNOVATIVE MEDICAL TECHNOLOGIES, INC
OMAHA, NE

866-560-1820
support@InnovativeMedTech.com
www.InnovativeMedTech.com